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<http://www.bijouna.jp/>

Bijouna

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Bijouna

Instruction Manual

Please read this manual before using the product.

Includes
warranty

Thank you for purchasing the Bijouna hot & cool facial massager. In order to learn how to use the product correctly, please read this instruction manual before use. After you have finished reading the manual, store it safely in a location where it can be referred to at any time.

Koizumi Seiki Corp.

Part names

KBE-1610

Facial massager

Cotton pad cover

Sets a cotton pad on the massager.

Head

(Titanium coating)

* Please remove the protective sticker before use.

Vents

Ground

Ground (sides)

Please use the product by holding (touching) this part with your bare hand.

COOL button

Display

HOT button

Power/selection button


(Hot & cool/cleansing/moisturizing care)

Operating/charging lamp

Back

Socket

AC adapter (KTH-9003)

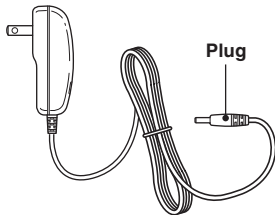Accessories

Moisture serum exclusively for this product
(sample size, 10 mℓ)

* Our KBE-C006 moisture serum (30 mℓ) exclusively
for this product is sold separately.

Charging the product


This product is for corded and cordless use.

You can use the product when it is connected to the AC adapter.

When you first use the product after purchase and when you have not used the product for a long time, charge it for about three hours.

1 Push the plug into the socket.

2 Plug the AC adapter into an outlet.
(When charging, the operating/charging lamp is illuminated in red.)

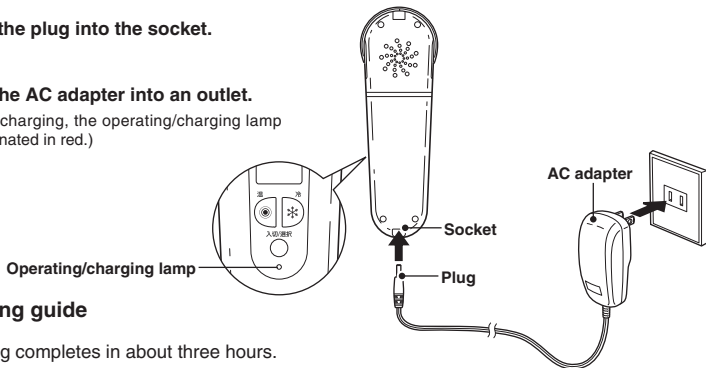Charging guide

- Charging completes in about three hours.
(The operating/charging lamp will turn off.)

Attention

- The product and AC adapter may feel warm when charging or during use. This is normal.
- Do not charge the product with any AC adapter other than the one exclusively for this product.

After charging ends

Unplug the AC adapter from the outlet and unplug the plug from the socket.

Before use

Warning

Stop using the product and consult with a medical specialist when you feel pain or you feel discomfort when using the product, or when your skin feels unusual after using the product. Otherwise, there is a risk of an accident, health problems, or skin problems.

Caution

Do not use this product if you have a skin injury or inflammation. Otherwise, there is a risk of skin problems or pain.

Remove any accessories such as earrings, piercings, and necklaces, and remove your contact lenses before using the product. Otherwise, there is a risk of an accident or other problems.

Modes

* Always observe the usage frequency and time. * When the product is powered, an operation confirmation beep sounds every five seconds.

HOT

* No longer than five minutes per treatment
Slowly warms your skin.

HOT cleansing

* No longer than five minutes per treatment, once a day
Removes microscopic dirt that cannot be cleaned with normal face washing.

COOL

* No longer than 90 seconds per treatment
Cools your skin to tighten your pores.

HOT moisturizing care

* No longer than five minutes per treatment, two or three times a week
Allows the beauty components of the moisture serum to penetrate into the outermost layer of the skin.

Suggested courses

Full course (two or three times a week)

step
1

HOT cleansing

step
2

HOT moisturizing care

step
3

COOL

Preparation before treatment

- Clean your skin by washing your face and removing any makeup, oil, and dirt.
- Remove any accessories and metallic objects on your body.

step
1

HOT

Preparation before treatment

step
2

COOL

- Remove any accessories and metallic objects on your body.

Hot & cool care course (daily)

HOT (no longer than five minutes)

Time per treatment

Auto power off - Approx. five minutes

Attention


Use the product in an ambient temperature of 10 to 40°C.

If the product is used in high temperatures, the head will get too hot and damage your skin.

1 Turn the power on and select "HOT" mode.

Press and hold the "Power/selection button" to turn on the power.
("HOT" mode is selected when the power is turned on.)

- The display changes as shown below until the temperature has stabilized.
When the display changes to the [Stabilized] state, you may start using the product.

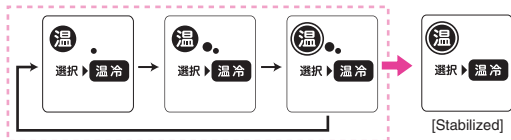2 Start the treatment.

- Hold the product so that you do not block its vents.
- Lightly press the head against your skin and move it slowly across your face by patting the skin.


* When the product is powered, an operation confirmation beep sounds every five seconds.

Attention

- Do not use the product continuously for over five minutes per treatment. Otherwise, there is a risk of skin problems.
- Do not slide the head against your skin.
- Do not stop the head on a single portion of your skin, move it slowly across your face by patting the skin.

- This product makes a buzzing noise. This is normal.

* Hold the product so that you do not block its vents.

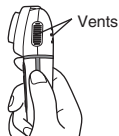- Forehead: From bottom to top
- Nose bridge: From bottom to top
- Cheeks: From center to outside
- Sides of nostrils: From top to bottom
- Underneath nose: From top to bottom
- Chin: From top to bottom

- * The "HOT" treatment is no more than five minutes per treatment. When five minutes has elapsed, the power automatically turns off.


After use

3 Press and hold the "Power/selection button" to turn off the power.

Attention

- **Do not continuously use the product.**

This places a strain on the skin and may cause skin problems or low-temperature burns.

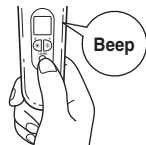HOT cleansing

Frequency


Time per treatment

(once a day/no longer than five minutes)

Auto power off
- Approx. five minutes

1 Place a cotton pad on the head.

- Place a cotton pad on the head and fasten it securely with the cotton pad cover.
- The head is large, so you will need large cotton pads.

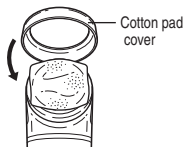2 Moisten the cotton pad with moisture water.

- For moisture water, use our KBE-C007 moisture water (120 mℓ) exclusively for this product (sold separately).
- Moisten the cotton pad to a degree where the moisturizer does not drip from it.
(So that the cotton pad becomes slightly transparent and the head can be seen)


3 Turn the power on and select "HOT cleansing" mode.

- ① Press and hold the "Power/selection button" to turn on the power.
("HOT" mode is selected when the power is turned on.)

- ② Press the "Power/selection button" once to select "HOT cleansing" mode.

- The display changes as shown below until the temperature has stabilized.
When the display changes to the [Stabilized] state, you may start using the product.

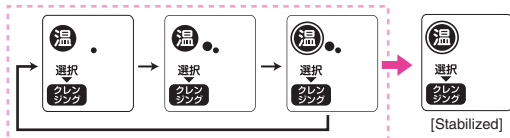4 Start the treatment.

① Hold the product so that your bare hand touches the ground.

Attention

- If you do not touch the ground with your bare hand during use, you cannot perform a "HOT cleansing" treatment. Always hold the product so that your bare hand touches either the left or right ground.

② Lightly press the head against your skin and slowly move from the center of your face to the outside.

* When the product is powered, an operation confirmation beep sounds every five seconds.

Attention

- Do not use the product continuously for over five minutes per treatment, once a day. Otherwise, there is a risk of skin problems.
- Do not stop the head on a single portion of your skin, use the product by moving the head around. Otherwise, there is a risk of skin problems or health problems.
- If the cotton pad becomes dry or it is hard to slide, add moisture water.
- This product makes a buzzing noise. This is normal.

After use

5 ① Press and hold the "Power/selection button" to turn off the power.
② Open the cotton pad cover and remove the cotton pad.


- Remove the cotton pad and wipe off any moisture water on the product including the head and cotton pad cover using a soft cloth or tissue.

Attention

- Do not continuously use the product.
This places a strain on the skin and may cause skin problems or low-temperature burns.
- Do not reuse any cotton pads that you have already used.
Otherwise, there is a risk of skin problems.

* If you will not perform the "HOT moisturizing care" treatment after you finish the "HOT cleansing" treatment, perform your normal facial routine such as applying moisture water and moisture serum to your skin.

* If you will continue on to the "HOT moisturizing care" treatment, first treat your skin with our KBE-C007 moisture water (120 ml) exclusively for this product (sold separately).

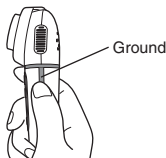- ① Forehead: From bottom to top ④ Sides of nostrils: From top to bottom
② Nose bridge: From bottom to top ⑤ Undereath nose: From top to bottom
③ Cheeks: From center to outside ⑥ Chin: From top to bottom

- The "HOT cleansing" treatment is no longer than five minutes once a day. When five minutes has elapsed, the power automatically turns off.

HOT moisturizing care

Frequency

Time per treatment

Auto power off
- Approx. five minutes

(two or three times a week/no longer than five minutes)

1 Evenly apply moisture serum to your entire face.




- For the moisture serum, use our KBE-C006 moisture serum (30 mℓ) exclusively for this product (sold separately).
- * This product includes our moisture serum (sample size, 10 mℓ) exclusively for this product.

2 Turn the power on and select "HOT moisturizing care" mode.

- ① Press and hold the "Power/selection button" to turn on the power.
("HOT" mode is selected when the power is turned on.)

- ② Press the "Power/selection button" twice to select "HOT moisturizing care" mode.

- The display changes as shown below until the temperature has stabilized.
When the display changes to the [Stabilized] state, you may start using the product.

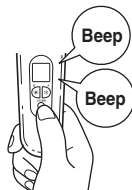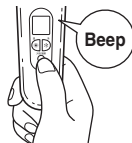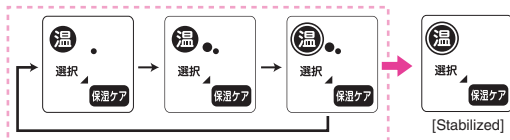
4 Start the treatment.

① Hold the product so that your bare hand touches the ground.

Attention

- If you do not touch the ground with your bare hand during use, you cannot perform a "HOT moisturizing care" treatment. Always hold the product so that your bare hand touches either the left or right ground.

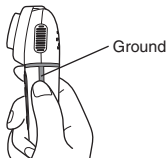② Lightly press the head against your skin and slowly move from the center of your face to the outside.

* When the product is powered, an operation confirmation beep sounds every five seconds.

Attention

- Do not use the product continuously for over five minutes per treatment, two or three times a week. Otherwise, there is a risk of skin problems.
- Do not stop the head on a single portion of your skin, use the product by moving the head around. Otherwise, there is a risk of skin problems or health problems.
- If the head is difficult to slide, add moisture serum.
- This product makes a buzzing noise. This is normal.

① Forehead: From bottom to top ④ Sides of nostrils: From top to bottom
② Nose bridge: From bottom to top ⑤ Underneath nose: From top to bottom
③ Cheeks: From center to outside ⑥ Chin: From top to bottom

- The "HOT moisturizing care" treatment is no more than five minutes per treatment, two or three times a week. When five minutes has elapsed, the power automatically turns off.

After use

5 Press and hold the "Power/selection button" to turn off the power.

- Wipe off any moisture serum on the product including the head and cotton pad cover using a soft cloth or tissue.

Attention

- Do not continuously use the product. This places a strain on the skin and may cause skin problems or low-temperature burns.

* If you will not perform any other treatments after you finish the "HOT moisturizing care" treatment, perform your normal facial routine such as applying moisture water and moisture serum to your skin.

COOL (no longer than 90 seconds)

Time per treatment

Auto power off - Approx. 90 seconds


Attention

Use the product in an ambient temperature of 10 to 40°C.

If the product is used in low temperatures, the head will get too cold and damage your skin.



1 Turn the power on and select "COOL" mode.

- ① Press and hold the "Power/selection button" to turn on the power.
("HOT" mode is selected when the power is turned on.)

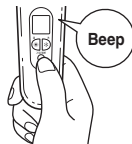- ② Press the "COOL button" to select "COOL" mode.

- The display changes as shown below until the temperature has stabilized.
When the display changes to the [Stabilized] state, you may start using the product.

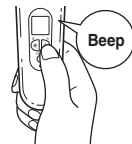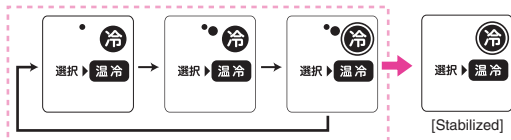2 Start the treatment.

① Hold the product so that you do not block its vents.

② Lightly press the head against your skin and move it slowly across your face by patting the skin.


* When the product is powered, an operation confirmation beep sounds every five seconds.

Attention

- Do not use the product continuously for over 90 seconds per treatment. Otherwise, there is a risk of skin problems.
- Do not slide the head against your skin.
- Do not stop the head on a single portion of your skin, move it slowly across your face by patting the skin.

• This product makes a buzzing noise. This is normal.

* Hold the product so that you do not block its vents.

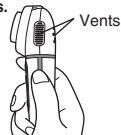- ① Forehead: From bottom to top ④ Sides of nostrils: From top to bottom
② Nose bridge: From bottom to top ⑤ Underneath nose: From top to bottom
③ Cheeks: From center to outside ⑥ Chin: From top to bottom


• The "COOL" treatment is no more than 90 seconds per treatment. When 90 seconds has elapsed, the power automatically turns off.

After use

3 Press and hold the "Power/selection button" to turn off the power.

Attention

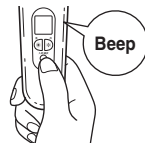- Do not continuously use the product.
This places a strain on the skin and may cause skin problems or low-temperature burns.

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